



# **User Manual**

J.Adamczewski-Musch, S.Linev, H.G.Essel GSI Darmstadt, Experiment Electronics Department

Produced: February 13, 2009, Revisions:

Titel: DABC: User Manual

Document	Date	Editor	Revision	Comment
DABC-user	2009-01-05	Hans G.Essel	1.0.0	First scetch

# **Contents**

2 DABC User Manual: Overview	5
2.1 User manual overview	. 5
2.1.1 Download and installation	. 5
2.1.2 Application working directory	. 5
2.1.3 Operation	. 5
3 DABC User Manual: Setup	7
3.1 Installing DABC	. 7
3.2 Setting up DABC user workspace	. 7
3.2.1 DABC data simulator	. 7
3.2.2 DABC MBS event server	. 7
3.2.3 PCI connected front-ends	. 7
4 DABC User Manual: GUI	9
4.1 GUI Guide lines	. 9
4.2 GUI Panels	. 10
4.2.1 Main DABC GUI buttons	. 10
4.2.2 DABC control panel	. 12
4.2.2.1 DABC controller buttons	. 13
4.2.3 Action in progress	. 13
4.2.4 MBS control panel	. 13
4.2.5 Combined DABC and MBS control panel	. 13
4.2.6 Command panel	. 14
4.2.7 Parameter table	. 14
4.2.7.1 Parameter selection	. 15
4.2.8 Monitoring panels	. 15
4.2.8.1 States	. 15
4 2 8 2. Rate meters	16

4.2.8.3 Histograms	16
4.2.8.4 Information	16
4.2.8.5 Logging window	16
4.3 GUI save/restore setups	18
5 DABC User Manual: MBS GUI	19
5.1 MBS event building	19
5.1.1 MBS setup	19
5.1.2 MBS control panel	19
5.1.2.1 MBS controller buttons	20
5.1.3 MBS command panel	21
5.2 MBS DIM parameters	21
5.2.1 MBS states	21
5.2.2 MBS rates	21
5.2.3 MBS infos	23
5.2.4 MBS tasks	23
5.2.5 MBS text	23
5.2.6 MBS numbers	23
5.3 Working directories	23
5.3.1 MBS configuration of DIM	23
6 DABC User Manual: DABC Application MBS	25
6.1 MBS event building with DABC	25
6.1.1 MBS setup	25
6.1.2 DABC setup	25
6.1.3 Combined DABC and MBS control panel	26
6.1.3.1 Combined DABC and MBS controller buttons	26
References	29
Index	29
	-

# DABC User Manual: Overview

[user/user-overview.tex]

#### 2.1 User manual overview

In the following sections we describe the usage on *DABC* applications.

#### 2.1.1 Download and installation

This starts with the installation described in section 3.1, page 7.

#### 2.1.2 Application working directory

We assume that all the necessary components of the application are ready. All plug-ins are implemented. Then we need mainly the setup files. These are typically in a working directory. More in section 3.2, page 7.

#### 2.1.3 Operation

Once we have the working directory setup we can begin. In the GUI part we describe the general functionality of the GUI which is common for most applications.

More specific functionality is described in the sections of the applications (as far as delivered with DABC).

As examples we take sometimes the *DABC* simulator plug-in which runs out of the box.

External application specific GUI add-ons cannot be described here.

# DABC User Manual: Setup

[user/user-setup.tex]

#### 3.1 Installing *DABC*

#### 3.2 Setting up *DABC* user workspace

To operate a *DABC* application one should create a dedicated working directory to keep all relevant files:

- Setup files for *DABC* (XML).
- Log files (text).

The GUI may run on a different directory with no access to the *DABC* working directory. There we place:

- Data files for startup panels (XML).
- Configuration files for GUI (XML).

Of cause, one can use one directory for all. The setup files are application specific. Therefore they are described in the application sections.

#### 3.2.1 DABC data simulator

#### 3.2.2 DABC MBS event server

5.1, page 19

#### 3.2.3 PCI connected front-ends

# DABC User Manual: GUI

[user/user-gui.tex]

#### 4.1 GUI Guide lines

The current *DABC* GUI is written in Java using the DIM software as communication layer. The standard part of the GUI described here may be extended by application specific parts. How to add such extensions is described in the programmer's manual. Typically they are started as prompter panels via buttons in the main GUI menu.

The standard part builds a set of panels (windows) according the parameters the DIM servers offer. Only services from one single DIM name server (node name specified as shell variable DIM\_DNS\_NODE) defining a name space can be processed.

The GUI needs no file access to the *DABC* working directory. However, user must have ssh (or rsh) access to the *DABC* (or *MBS*) master node. Currently the GUI must run under the same account as the *DABC*. In spectator mode (no commands) the GUI may run under different account. Master node must have remote access to all worker nodes. The user's ssh settings must enable remote access without prompts.

The layout of the GUI can be adjusted to individual needs. It is strongly recommended to save these settings to see the same layout after a restart of the GUI. The GUI can be restarted any time. *DABC* and *MBS* systems continue without GUI.

The GUI is packed in a JAR file xgui.jar which can be moved everywhere. The file must be in the CLASSPATH paths together with the pimdim/classes path. The GUI is invoked by java-Xmx200m xgui.xGui. For standard installations short cuts are provided like:


**dabc**: Start standard *DABC* GUI with *DABC* control panel.

**mbs**: Start GUI with MBS control panel only. This is for controlling a stand-alone running MBS.

**dabs**: Start GUI with combined MBS and DABC control panel. This is for using DABC as event builder with MBS front-ends.

moni : Start GUI without active elements. Works then as monitor.

guru: Starts GUI with all buttons for all functions.

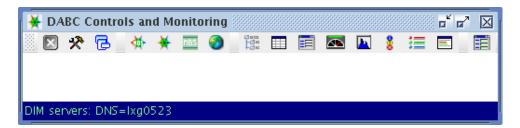

**Figure 4.1:** Main toolbar buttons.

#### 4.2 GUI Panels

Fig. 4.1, page 10 shows the main menu of *DABC* (minimal view). The GUI as it comes up is divided in three major parts: one sees on top a toolbar with icon buttons. Most of these open other windows. The dark line at the bottom shows a list of active DIM servers. The other windows are placed in the white middle pane. The functions of the buttons and the invoked panels is described in the next sections. Depending on the application some buttons may be not seen, additional ones may show up. If one does not work with *MBS* plug-ins the control panels for *MBS* are of cause not useful.

Fig. 4.2, page 11 shows a more typical view of a running *DABC*. In general, all panels (including the GUI itself) can be closed and reopened any time.

#### 4.2.1 Main DABC GUI buttons

Quit GUI. Will prompt (RET will quit). The *DABC* will continue to run. The GUI may be started anywhere again. In case you saved the layout (recommended, see 4.3, page 18) and you start the GUI from the same directory it will look pretty much the same as you left it.

- Test, shell script
- Save settings: window layout, record attributes, parameter selection filters. Details see 4.3, page 18. Note that the content of the control panels must be saved by similar buttons in these panels.
- Open *DABC MBS* control panel, see 6.1.3, page 26.
- Open DABC control panel, see 4.2.2, page 12.
- Open MBS control panel, see 5.1.2, page 19.
- Refresh. Update DIM service list from DIM name server. New parameters and commands are added. Existing ones are not removed, but only deactivated. Parameters and Commands are sorted alphabetically by name. All panels are updated. This update function is called automatically whenever the frontend systems might have added or deleted DIM services. In normal operation there is no need to update manually.
- Open command panel (4.2.6, page 14).
- Open parameter table (4.2.7, page 14).
- Open parameter selection panel (4.2.7.1, page 15).
- Open rate meter panel (4.2.8, page 15).

4.2. GUI Panels


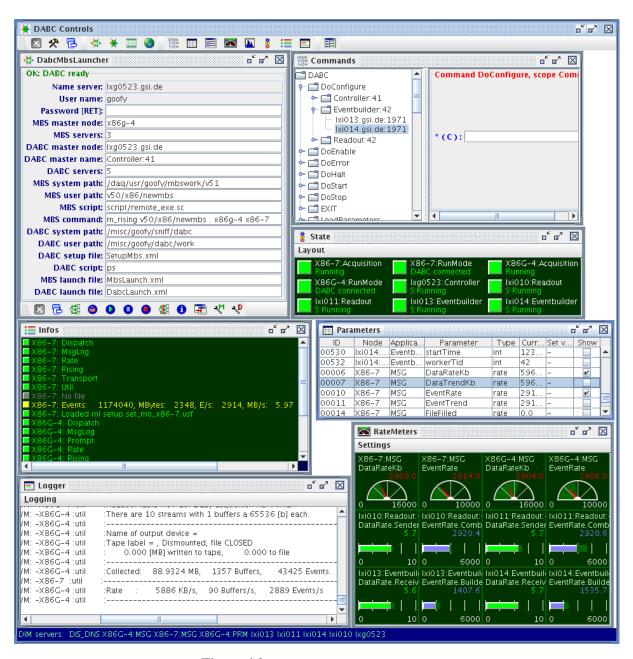


Figure 4.2: More typical full screen view.

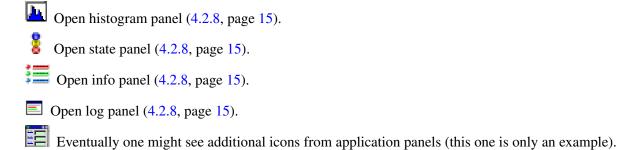

The three control panels (*DABC*, *MBS*, combined *DABC* and *MBS*) are used depending on the application to be controlled. Eventually an application provides additional specific control panels.

#### 4.2.2 DABC control panel

The standard *DABC* control panel is shown in 4.3, page 12. As mentioned already some applications may provide their own control panels like the *MBS* applications (see section 5.1.2, page 19). But most of the buttons are very common. From left to right they startup a system, configure it, start data taking, pause data taking, stop tasks, shut down. At the very left we see a save button, at the right a shell execution button. Values are read from file <code>DabcControl.xml</code> (default, may be saved/restored to/from other



**Figure 4.3:** DABC controller panel.

file, see 4.3, page 18).

```
<?xml version="1.0" encoding="utf-8"?>
<DabcLaunch>
<DabcMaster prompt="DABC Master" value="node.xxx.de" />
<DabcName prompt="DABC Name" value="Controller:41" />
<DabcUserPath prompt="DABC user path" value="myWorkDir" />
<DabcSystemPath prompt="DABC system path" value="/dabc" />
<DabcSetup prompt="DABC setup file" value="SetupDabc.xml" />
<DabcScript prompt="DABC Script" value="ps" />
<DabcServers prompt="%Number of needed DIM servers%" value="5" />
</DabcLaunch>
```

**DabcMaster:** Node where the master controller shall be started. Can be one of the worker nodes.

**DabcName:** A unique name inside *DABC* of the system.

**DabcUserPath:** User working directory. The GUI does not need to have access to the filesystem.

**DabcSystemPath:** Path where the *DABC* is installed.

DabcSetup: Setup file name.

**DabcScript:** Command to be executed in an ssh at the master node.

**DabcServers:** Number of workers and controllers. This information is minimum for the GUI to know when all *DABC* nodes are up. The GUI waits until this number of DIM servers is up and running.

The name server name is translated from shell environment variable DIM\_DNS\_NODE, the user name from shell environment variable USER. Password can be chosen when the first remote shell script is executed (which itself is protected by user password). All following commands then need this password.

4.2. GUI Panels

#### **4.2.2.1** *DABC* controller buttons

Save panel settings to the file Control file. If you choose a name different from the default you must set a shell variable to it to get the values from that file (see 4.3, page 18).

Startup all tasks. Executes a *DABC* script dabcstartup.sc via ssh on the master node under user name. Then it waits until the number of DIM servers expected are announced. A progress panel pops up during that time (see 4.2.3, page 13). When the servers are up the main GUI Update is triggered building all panels from scratch according the parameters offered by the servers.

Configure. Executes state transition command Configure on master node and waits for the transition. All plug-in components are created. Then execute Enable. Waits until all workers go into Ready state. Now the *DABC* is ready to run. Triggers the main GUI Update.

- Start acquisition. Executes Start command. All components go into running state Running.
- Pause acquisition. Executes Stop command. All components go into standby state Ready.
- Halt acquisition. Executes Halt command. This closes all plug-ins. States go into Halted. Next must be shut down or configure.
- Exit all processes by EXIT commands. After 2 seconds trigger the main GUI Update.
- Shut down all processes on all nodes by script. This is the hard shut down.
- ssh shell script execution on master node.

#### 4.2.3 Action in progress



Figure 4.4: Launching progress.

When starting up, configure or shut down the GUI has to wait until the front-ends have completed the action. During that time a progress window similar to the one shown in Fig. 4.4, page 13 pops up. Please wait until the popup disappears.

#### 4.2.4 MBS control panel

To control and monitor a stand-alone *MBS* system a dedicated control panel is provided by the *MBS* application. This panel is described in the *MBS* section 5.1.2, page 19.

#### 4.2.5 Combined *DABC* and *MBS* control panel

To control and monitor *MBS* front-ends with *DABC* event builders a dedicated control panel is provided by the *MBS* application. This panel is described in the *MBS* section 6.1.3, page 26.

#### 4.2.6 Command panel

The control system of *DABC* and/or the application specific plug-ins can define commands. These commands are encoded as DIM services including a full description of arguments. Therefore the GUI can build up at runtime a command tree and provide the proper forms for each command. Commands are executed in all components of *DABC*.

The DABC naming convention for commands and parameters defines four main name fields separated by slashes:

- 1. DIM server name space (example DABC)
- 2. Node (example dabcWorker1)
- 3. Application (example Combiner)
- 4. Name (example doConfigure)

Example: DABC/dabcWorker1/Combiner/doConfigure. Fig. 4.5, page 14 shows on the left side the command

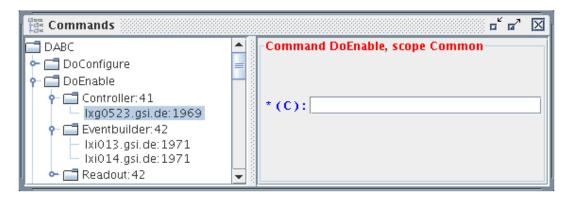


Figure 4.5: Command panel.

tree. The tree is built from name, application, nodes. Double click (or RETURN) on a treenode executes the command on all treenodes below. A click on a command opens at the right side the argument panel. Entering argument values and RETURN executes the command. In the example shown in the figure doble click on doEnable would execute that command on three nodes. Double click on Eventbuilder would execute only on two nodes.

#### 4.2.7 Parameter table

*DABC* parameters are DIM services as the commands. The naming convention is the same. The server providing parameters can be make them (no)visible and (un)changable. *DABC* defines some special parameter types having a data structure and a specific interpretation like a rate parameter having a value, limits, a color, and a graphic presentation. A rate parameter is assumed to be changed and updated regularly. The GUI displayes these special parameters in dedicated panels. Parameters are used in all components of *DABC*. The central place for all parameters

Ⅲ Parameters					X			
ID	Node	Application	Parameter	Туре	Current	Set value	Show	
00308	lxi010:1970	Readout: 42	CtrlPoolSize.BnetPlugin	int	2097152			
00309	lxi010:1970	Readout:42	DABCVersion	char	DABC C	_		
00310	lxi010:1970	Readout:42	DataRate.Sender	rate	0.0	_	V	
00311	lxi010:1970	Readout: 42	EventBuffer.BnetPlugin	int	524288			트
00312	lxi010:1970	Readout: 42	EventPoolSize.BnetPlugin	int	4194304			
00313	lxi010:1970	Readout: 42	EventRate.Combiner	rate	0.0	_	V	
00314	lxi010:1970	Readout: 42	InfoMessage	info	State ma	-	~	¥


**Figure 4.6:** Parameter table.

eters in the GUI is the parameter table as shwon in Fig. 4.6, page 14. The parameter table holds all parameters

4.2. GUI Panels

which are marked by the provider to be visible. The parameter values can be changed in the Set value column if no minus sign is there in which case the provider does not grant modification. The buttons in the Show column indicate if the parameter is shown in some graphics panel. It can be removed from or added to this panel by the buttons. The table can be ordered by columns (click on column header). The column width can be adjusted and is saved/restored by main save button (see 4.3, page 18).

#### 4.2.7.1 Parameter selection

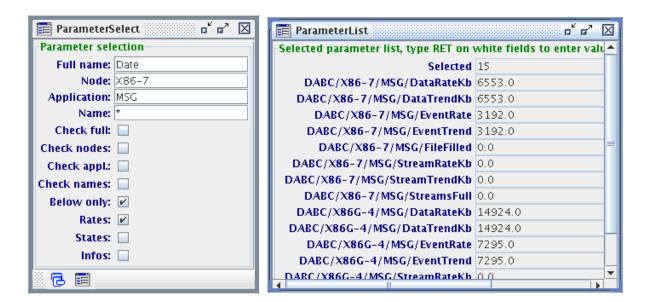

**Figure 4.7:** Parameter selection panel and selected parameter list.

To get a more selective view on the parameters one can specify filters in the panel shown at the left side of Fig. 4.7, page 15. Text substrings for each of the four name fields can be specified as well as a selection of record types. Values can be saved (see 4.3, page 18). With the check boxes the filter function for each of these can (de)activated. The parameter list at the right window in Fig. 4.7, page 15 shows only the parameters matching all filters.

If the data field is white the parameter can be changed. This cannot be done in place because the parameter might be updated in the mean time. Instead press RETURN in the field. A prompter will pop up to enter the value.

#### 4.2.8 Monitoring panels

As already mentioned the *DABC* provides definitions of special purpose DIM parameters. These Records can be recognized by the GUI and are handled in appropriate way. Currently there are

- o States
- o Rates
- o Histograms
- o Infos

#### 4.2.8.1 States

States are records having a number for severity (0 to 4), a color, and a brief state description (see Fig. 4.8, page 16). Of cause the states of the *DABC* state machine are shown as states. Application plug-ins may use this kind of records also for other information.



Figure 4.8: States.

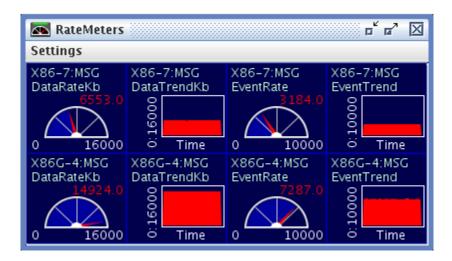


Figure 4.9: Rates.

#### **4.2.8.2** Rate meters

All rate meters are displayed in the meter panel, Fig.4.9, page 16. Meters can be removed in the parameter table (See Fig. 4.6, page 14) with the Show buttons like the other graphical parameters. Saving the setup, the visibility will be preserved.

On the left side in Fig. 4.10, page 17 the Settings menu is shown. It affects all items in the panel. One can Zoom (toggle between large and normal view), change the number of columns, change the display mode, toggle Autoscale, and set limits (applied to all meters).

Besides that each individual item can be adjusted by right mouse button. The context menu is shown on the right. All changes done individually are changing the defaults! The global changes can be overwritten by these defaults. All settings are saved with the setup and restored on GUI startup (see 4.3, page 18).

#### 4.2.8.3 Histograms

Histogram panels are handled in pretty much the same way as the rate meters. All histograms are displayed in the histogram panel, Fig.4.11, page 17. Histograms can have arbitrary size set in Layout menu.

#### 4.2.8.4 Information

Information records mainly display one line of text with a color (see Fig. 4.12, page 17).

#### 4.2.8.5 Logging window

Fig. 4.13, page 18 show the logging window.

4.2. GUI Panels

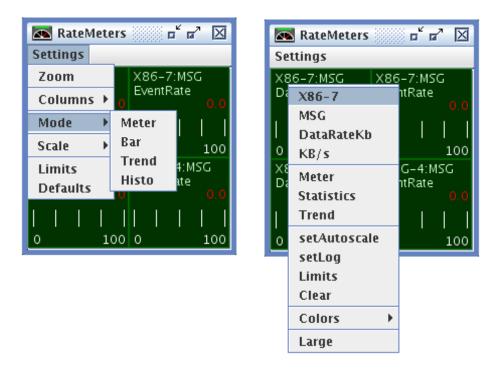


Figure 4.10: Steering menus.

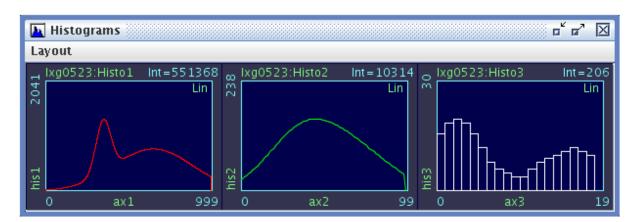


Figure 4.11: Histograms.

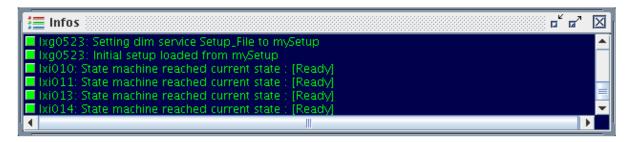


Figure 4.12: Info.

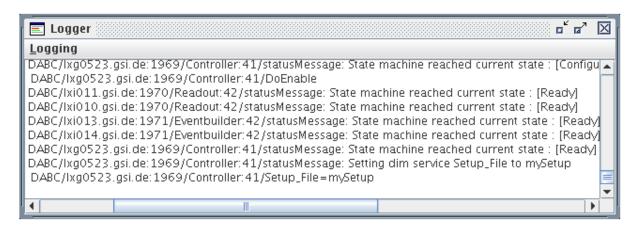


Figure 4.13: Logging.

#### 4.3 GUI save/restore setups

There are several setups which can be stored in XML files and are retrieved when the xGUI is started again. The file names can be specified by shell variables.

DABC\_CONTROL\_DABC: Values of *DABC* control panel. Saved by button in panel.

Default DabcControl.xml. Filename in panel itself.

DABC\_CONTROL\_MBS: Values of MBS control panel. Saved by button in panel.

Default MbsControl.xml. Filename in panel itself.

DABC RECORD ATTRIBUTES: Attributes of records. Saved by main save button.

Default Records.xml.

DABC PARAMETER FILTER: Values of parameter filter panel. Saved by main save button.

Default Selection.xml.

DABC GUI LAYOUT: Layout of frames. Saved by main save button.

Default Layout.xml.

## DABC User Manual: MBS GUI

[user/user-gui-mbs.tex]

#### 5.1 *MBS* event building

#### **5.1.1** *MBS* setup

Any MBS system can be controlled by the DABC GUI. It can run in two operation modes: with MBS event builder or DABC event builder (see 6.1, page 25. The first case means a standard MBS system.

To control a standard *MBS* nothing has to be done by the user on the *MBS* side. The node running the GUI must get granted rsh access at least to the *MBS* node where the prompter shall run. **Note**, however that in the user's *MBS* startup file (typically startup.scom) the m\_daq\_rate task must be started as last task (this is probably the case already). This task calculates the rates. The GUI waits for this task after execution of the startup file. Because *MBS* has no states there is no other way to know when the startup has finished. Of cause, the *MBS* itself must have been build with the DIM option (since version v5.1). Central log file is written as usual. Optionally one can provide a text file with specifications which parameters shall appear on screen (see 5.3.1, page 23).

For the standard MBS control one needs no DABC installation. The GUI jar file is sufficient. DIM must be installed.

#### **5.1.2** *MBS* control panel

Fig. 5.1, page 20 shows the panel to be used to control a standard MBS. The values are restored from file MbsControl.xml (default, may be saved to other file, see 4.3, page 18).

```
<?xml version="1.0" encoding="utf-8"?>
<MbsLaunch>
<MbsMaster prompt="MBS Master" value="node-xx" />
<MbsUserPath prompt="MBS User path" value="myMbsDir" />
<MbsSystemPath prompt="MBS system path" value="/mbs/v51" />
<MbsStartup prompt="MBS startup" value="startup.scom"/>
<MbsShutdown prompt="MBS shutdown" value="shutdown.scom"/>
<MbsCommand prompt="Script command" value="whatever command" />
<MbsServers prompt="%Number of needed DIM servers%" value="3" />
</MbsLaunch>
```

**MbsMaster**: Lynx node where the *MBS* prompter is started.

**MbsUserPath**: *MBS* user working directory. The GUI need not to have access to that filesystem. **MbsSystemPath**: Path on Lynx where the *MBS* is installed. GUI needs no access to this path.



Figure 5.1: MBS controller.

**MbsStartup**: The user specific *MBS* startup command procedure, typically startup.scom, located on user path.

MbsShutdown: The user specific \$MBS\$ shutdown command procedure, typically shutdown.scom, located on user path.

**MbsCommand**: With RET an *MBS* command in executed (on current node). The shell script button executes this string as rsh command on master node.

**MbsServers**: Number of nodes plus prompter. This information is minimum for the GUI to know when all *MBS* nodes are up. The GUI waits until this number of DIM servers is up and running.

That file can be created from within the GUI in the MBS controller panel. Enter all values necessary, and store them.

#### **5.1.2.1** *MBS* controller buttons

Save panel settings, see 4.3, page 18.

Execute script prmstartup.sc at master node. Starts prompter, dispatchers and message loggers and waits until they are up. Trigger the main Update. A progress panel pops up during that time (see 4.2.3, page 13).

Execute script dimstartup.sc at master node. Starts dispatcher and message logger for single node *MBS*. Trigger the main Update.

- Configure. Execute user's MBS startup procedure in prompter (dispatcher). Trigger the main Update.
- Start acquisition. Execute Start acquisition.
- Pause acquisition. Execute Stop acquisition.
- Halt acquisition. Execute user's *MBS* shutdown procedure in prompter. Prompter, dispatcher and message loggers should still be running.
- Shut down all. Execute script prmshutdown.sc at master node. After 2 seconds trigger the main Update.
- Show acquisition. Output in log panel.
- Shell script executes command on master node.

#### **5.1.3** *MBS* command panel

Commands of a 🖂				
← □ MbsTransport	•	Command OpenFile, scope MBS		
► 🗂 CloseFile		name (C, OPT):		
← 📹 ConnectRfio		tape (L, OPT):		
► 🗂 DisableTcp	3	sizemb (L, OPT):		
► 🗂 DisconnectRfio		number (L, OPT):		
→ 🛅 DismountTape		first (L, OPT):		
EnableDabc		inhead (C, OPT):		
EnableTcp	1	outhead (C, OPT):		
ritTape ☐ InitTape ☐ MountTape		-auto 🗌		
→ ☐ OpenFile		-prompt 🗌		
∽ ☐ SetFileheader		-edit □		
∽ 🗂 ShowFile	=	-disk		
∽ 🗂 ShowTape		-rfio		
⊶ 🗂 TypeEvent		-multi 🔲		
⊶ 🗂 TypeStream				
← 🗂 MbsUtil		-first 📗		
► 🗂 MbsVme_serv	<b>-</b>			

Figure 5.2: Command panel.

Fig. 5.2, page 21 shows on the left side the command tree. Double click (or RETURN) on a command executes the command. The top tree level is the executing *MBS* task, below that are the commands, and the master node (prompter node) is the only node below each command. However, command is sent to the prompter node, but executed on the current node which is displayed in the info panel (see Fig. 5.4, page 22). Click on a command opens at the right side the argument panel. Entering argument values and RETURN executes the command. Only the *MBS* commands of the running tasks are shown. Fig. 5.3, page 22 shows that only dispatcher and prompter are up and therefore only their commands are seen. Fig. 5.4, page 22 shows in addition the commands of util and transport after configuration.

### 5.2 MBS DIM parameters

#### **5.2.1** *MBS* states

Acquisition/State Running | Stopped BuildingMode/State Delayed | Immediate EventBuilding/State Working | Suspended FileOpen/State File open | File closed

RunMode/State DABC connected I MBS to DABC | Transport client I MBS standalone

**SpillOn/State** Spill ON | Spill OFF **Trigger/State** Master | Slave

#### **5.2.2** *MBS* rates

MSG/DataTrendKb KBytes/s as trend MSG/EventRate Events/s MSG/EventTrend Events/s as trend

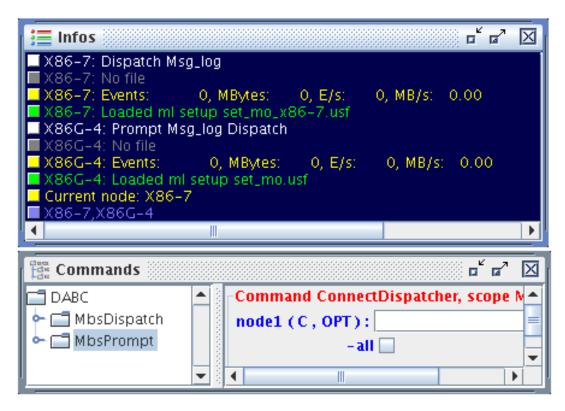


Figure 5.3: Info and command panel.

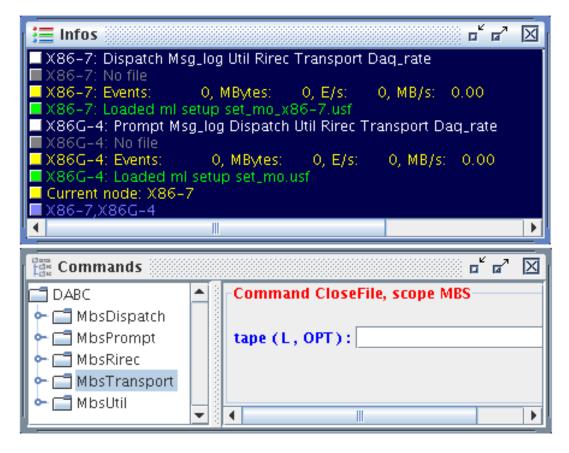


Figure 5.4: Info and command panel.

MSG/StreamRateKb Stream server Kbyte/s
MSG/StreamTrendKb Stream server Kbyte/s as trend
MSG/FileFilled File filled in percent
MSG/StreamsFull Number of full streams in percent

#### **5.2.3** *MBS* infos

Shown in info window.

MSG/eFile Name of file.

MSG/ePerform Events, MBytes, Events/s and MBytes/s.

**MSG/eSetup** Name of setup file loaded.

PRM/Current Current command execution node (master node only).

**PRM/NodeList** List of nodes (master node only).

#### **5.2.4** *MBS* tasks

Task list is shown in info window (name slightly different).

Dispatch MsgLog ReadMeb Collector Transport EventServ Util ReadCam EsoneServ StreamServ Histogram Prompt Rate SMI Sender Receiver AsynchReceiver Rising TimeOrder VmeServ

#### **5.2.5** *MBS* text

MSG/GuiNode Node where GUI runs
MSG/Date Date as written in file header
MSG/Run Run ID as written in file header
MSG/Experiment Experiment as written in file header
MSG/User Lynx user name as written in file header
MSG/Platform CPU platform

#### 5.2.6 MBS numbers

MSG/BufferSize
MSG/Buffers collected so far.
MSG/Events collected so far.
MSG/FileMbytes written in file.
MSG/FlushTime
MSG/MBytes collected so far.
MSG/StreamKeep
MSG/StreamMbytes
MSG/StreamScale
MSG/StreamSync

#### 5.3 Working directories

#### 5.3.1 MBS configuration of DIM

Optional file .guirc in the MBS working directory specifies which rate meters and states shall appear in the GUI. Upper limits of the rate meters can be specified. This file can be copied from \$MBSROOT/set. Only the parameters which are in this file are optional.

## This file controls the rate meter and state appearance.

```
## File name must be .guirc and in the MBS working directory.
## The value numbers are the maximum values for rate meters
## Colons only if value is specified!
## Node names must be uppercase, * wildcards all
## RunMode shows if MBS is set into DABC event builder mode
## and client is connected.
##====== All nodes:
##---- Rates:
               : 10000.
* EventRate
              : 10000.
* EventTrend
              : 16000.
* DataRateKb
* DataTrendKb
               : 16000.
* StreamRateKb : 16000.
* StreamTrendKb : 16000.
# ++ File filling status in percent, typically only on one node (transport)
#* FileFilled : 100.
* StreamsFull : 100.
##---- States:
# ++ Delayed or immediate event building:
* BuildingMode
# ++ Current eventbuilding running or suspended:
* EventBuilding
# ++ Running mode, stand alone or DABC:
* RunMode
# ++ Shows spill signal:
#* SpillOn
# ++ Shows if file is open, typically only on one node (transport)
#* FileOpen
#* Trigger
##===== Node XXX
#XXX EventRate : 10000.
#XXX DataRateKb : 16000.
#XXX RunMode
#XXX FileOpen
#XXX FileFilled : 100.
#XXX SpillOn
#XXX EventTrend : 10000.
#XXX DataTrendKb : 16000.
```

# **DABC** User Manual: **DABC** Application **MBS**

[user/user-app-mbs.tex]

#### 6.1 *MBS* event building with *DABC*

To run *MBS* front-ends with *DABC* nodes as event builders some modifications of the *MBS* setup files must be done. For the *DABC* side setup files must be provided.

#### **6.1.1** *MBS* setup

When we want to use *DABC* nodes as event builders, we need a different setup on the *MBS* side. We assume that we have more than one *MBS* node. Such multi-node system is controlled by an *MBS* prompter running on one node. The setup has to be changed such that all nodes run as if they are stand alone. That means that each node runs the Readout - Collector - Transport chain. The *DABC* event builders connect to these transports. The *MBS* buffer size should be set to the stream size and the number of buffers per stream must be set to one.

#### **6.1.2** *DABC* setup

On the DABC user working directory we need configuration files.

Summary of parameters:

**MbsFileName** File name for list mode data file (LMD). Overwritten by command.

MbsFileSizeLimit File closes when size is reached, and new file opens.

BufferSize Should match MBS buffer size.

MbsServerKind Transport | Stream.

**MbsServerPort** Port number transport (6000).

**MbsServerName** *MBS* node of transport.

NumInputs Number of MBS channels for one combiner.

**DoFile** Provide output file. **DoServer** Provide server.

These parameters are used to configure an event generator:

NumSubevents FirstProcId SubeventSize Go4Random The following example configuration file \$DABCSYS/applications/mbs/Combiner.xml shows how to configure one combiner module reading from three *MBS* transport servers. Any running program is described by a Context, in this case named Combiner. The Run specifies the library, start function and log file name. We have one module Combiner (could be different name from context) with three input ports and two output ports.

Now we can use the combined controller panel to startup MBS and DABC.

#### 6.1.3 Combined DABC and MBS control panel

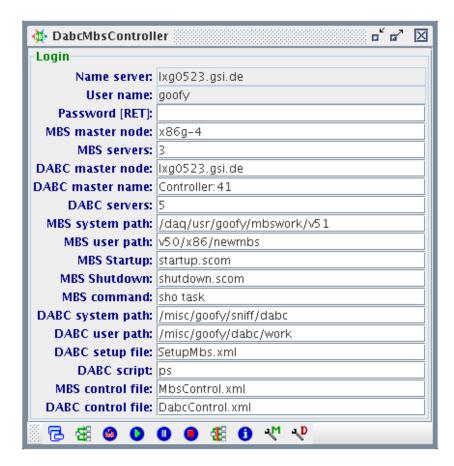


Figure 6.1: Combined DABC and MBS controller.

This panel shown in Fig. 6.1, page 26 is simply a superposition of the single ones.

#### 6.1.3.1 Combined DABC and MBS controller buttons

Save panel settings, see 4.3, page 18.

Execute script dabcstartup.sc at *DABC* master node. Starts DIM servers. Execute script prmstartup.sc at *MBS* master node. Starts prompter, dispatchers and message loggers. Waits for all components. A progress panel pops up during that time (see 4.2.3, page 13). If all components are up trigger the main Update.

Configure. Execute user's *MBS* startup procedure in prompter. Executes state transition command Configure on *DABC* master node and wait for the transition. All plug-in components are created. Then execute Enable. If all components are up trigger the main Update.

Start MBS acquisition, then executes DABC Start command. All components go into running state Running.

Pause acquisition. Execute *MBS* stop acquisition. Execute *DABC* Stop command. All components go into standby state Ready.

Halt acquisition. Executes *DABC* Halt command. This closes all plug-ins. States go into Halted. Execute user's *MBS* shutdown procedure in prompter. Prompter, dispatcher and message loggers should still be running. Next must be shut down or configure. After two seconds trigger the main Update.

Shut down all. Execute EXIT command on all *DABC* nodes. Execute script prmshutdown.sc at *MBS* master node. After two seconds trigger the main Update.



Shell script for *MBS* master node.

Shell script for *DABC* master node.

# **Index**

```
TODO
dabcsetupfiles, 26
Install, 7
Setup, 7
```